


520

Assisted Recline  
Tilt Chair featuring  
**contemporary styling**

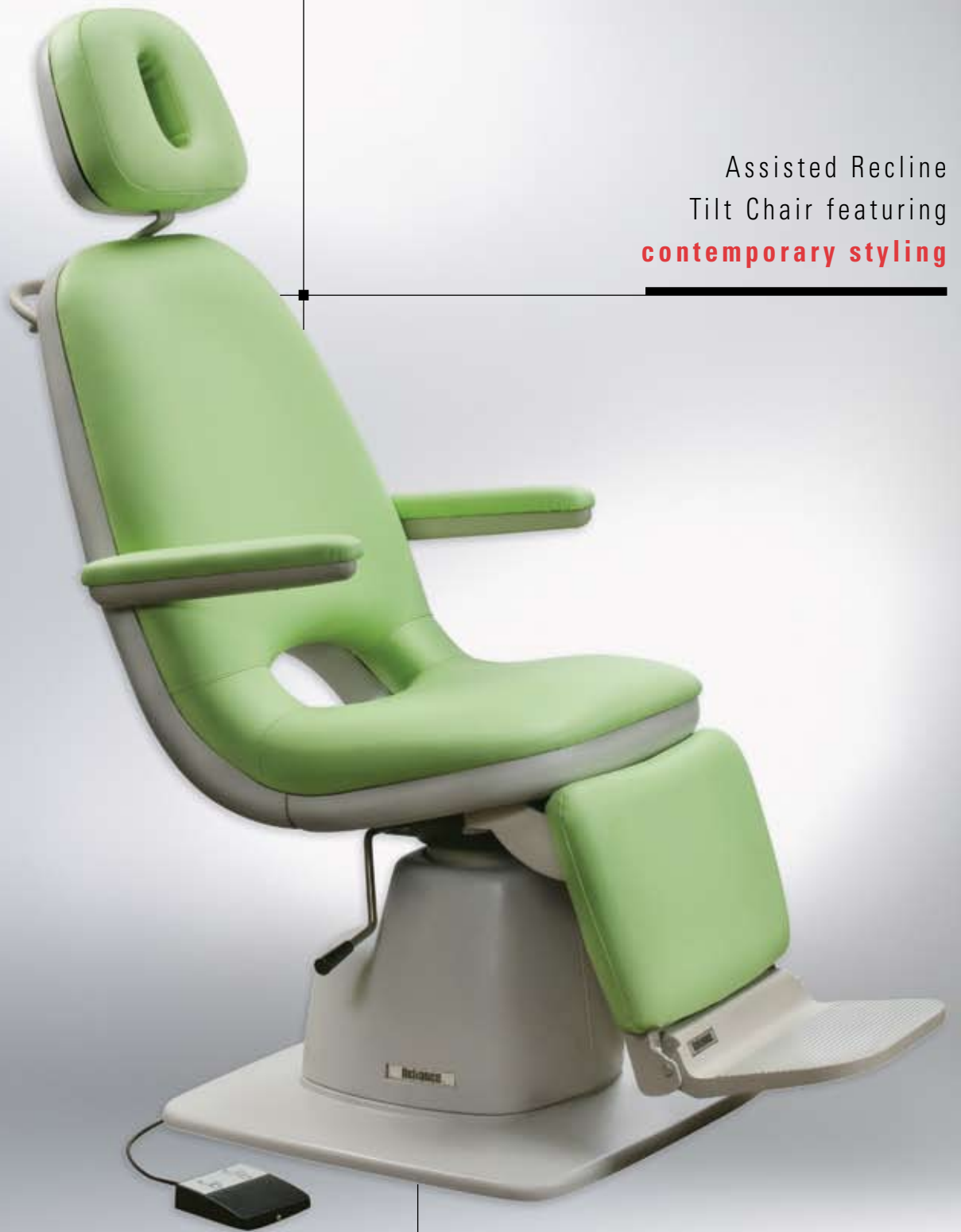

**Reliance**<sup>®</sup>

# 520

## Manual Assisted Recline



### Manual Reclining Top

A recline bar, which is easily grasped from both sides of the chair, controls the assisted recline of the chair. The chair reclines to the maximum angle of 40 degrees.

### Power Hydraulic Base

A corded footswitch controls the raise and lower functions of the base allowing you to maintain contact with the patient and eliminate wasted movement.

### Dual Rotation Locks

Easily accessible levers located on both sides of seat enable practitioner to rotate chair 360 degrees without changing position.

### Large Hinged Footstep

The 6"x14" footstep can support patient's full weight and offers safe ingress and egress. Footstep tilts out of the way when desired.

### Headrest

Large, padded headrest is fixed in place to keep the patient comfortable in the reclined position.

### Armrests

Upholstered arms rotate upward for convenient ingress and egress from either side of the chair. The armrests are upholstered and padded for patient comfort.

### Finishes and Upholstery

Hand-sewn upholstery with high density foam holds its shape for years. Cushions remove for cleaning and replacement. Match any office décor with our selection of 14 easy-to-maintain vinyl colors. In addition, you may choose one of 165 custom colors. Base is available in #20 Pearl. All materials are fire retardant. To request a color selection card, contact your dealer or call Reliance at 1-800-735-0357.

### Dimensions and Specifications

Low base height only: 21 1/2" lowered, 32 1/2" raised. Base: 22 1/8" x 24 1/4". Width: 28". Depth upright: 43 1/2". Reclined: 59 1/2". Seat: 21 1/2" (wide) x 15" (deep). Distance between armrests: 21 5/8". Approximate shipping weight: 310 lbs. Electrical: (Standard) 120 volts, 60 Hz, 3.0 amperes; (Optional) 230 volts, 50 Hz, 3.0 amperes. CSA Approved.



Precision by Tradition

Reliance® Medical Products • 3535 Kings Mills Road • Mason, Ohio 45040-2303 • 1-800-735-0357 • [www.reliance-medical.com](http://www.reliance-medical.com)

© 2007 Reliance Medical Products, Inc. All rights reserved. Printed in USA with soy inks on recycled paper. 9/07